

DEFIGARD TOUCH 7 PHYSIOGARD TOUCH 7

Defibrillator / Monitor

Maximum performance, minimum weight, touch screen control: the power to save lives



SCHILLER
The Art of Diagnostics

DEFIGARD TOUCH 7 PHYSIOGARD TOUCH 7

A compact combination of state-of-the-art defibrillation technology and comprehensive monitoring functions

The DEFIGARD Touch 7 and PHYSIOGARD Touch 7 (DGT7/PGT7) are indispensable tools for rescuers. In a defibrillator/monitor or only monitor version, the device is extremely compact and offers the latest defibrillation technology in combination with comprehensive monitoring functions. It is the first emergency monitor/defibrillator equipped with a touch screen – making it the most intuitive device on the market – as well as with the latest data transmission technology.



**DEFIGARD TOUCH 7
MONITOR / DEFIBRILLATOR**



**PHYSIOGARD TOUCH 7
MONITOR**

UNPARALLELED TOUCHSCREEN

The DGT7/PGT7 is the first emergency monitor/defibrillator with a touch screen. Thanks to the new software environment Touch'n Save, this device offers unmatched user-friendliness. Moreover, the exceptional display quality enables optimal display readability at all times. The large 7-inch touchscreen is particularly bright, and its high resolution guarantees unparalleled user-friendliness. Thanks to Projected Capacitive Touch technology, the display enables highly responsive operation, even when wearing gloves.

CAN BE USED ANYWHERE

Thanks to SCHILLER's extensive experience in the rescue sector, the DGT7/PGT7 has been developed so that it can also be used in extreme conditions. It is the ideal companion for all land and water interventions.

Through its lightweight structure, it also offers new possibilities for air rescue operations. Its low weight – three to four times below that of competing products – ensures significant fuel saving when used on an air plane or a helicopter.

With its mounting bracket – which also serves as battery charger – the DGT7/PGT7 is resistant to shocks and vibrations, thus complying with land and air rescue operation standards.

CONNECTED AT ALL TIMES

The device can perform high speed data transmission at all times. Depending on the location, transmission can be performed using Wi-Fi or 3G.

The safety standards WPA2 Enterprise guarantee secure data protection.

PRACTICAL AND SOPHISTICATED

In the functional carrier bag, each accessory has its place, is easily accessible and can readily be used, as each part is already connected. Moreover, the device can be tilted for optimal convenience.

OPERATIONAL AT ALL TIMES

The DGT7/PGT7 is ready to use at all times, thanks to regularly performed self-tests and to a RFID transponder that ensures the monitoring of the defibrillator electrodes' expiry date.



**ONLY 2.4 KG.
LIGHTWEIGHT, COMPACT
AND USER-FRIENDLY:
YOUR PERFECT COMPANION!**

ADAPTIVE USER INTERFACE



PHYSICIAN MODE This display, showing values and waveforms as well as a manual defibrillation function, gives the physician a comprehensive overview of the patient.



RESCUE MODE With a simple swipe gesture, the user can move from the physician display to the rescue display. The values are displayed numerically and are well visible, even from a distance. Moreover, a semi-automatic defibrillation function is also available for all users.

ECG

The DGT7 / PGT7 enables continuous ECG monitoring as well as 12-lead ECG diagnostic recording. ECG electrodes are conditioned so as to be used in emergency situations.

CAPNOGRAPHY

The EtCO₂ value is a parameter to check blood circulation and deduct the efficiency of CPR.

The DGT7 / PGT7 offers Mainstream or Sidestream measurement using MASIMO sensors to display EtCO₂, Respiration Rate and Apnea information.

SPO₂ , SPCO AND SPMET

The Masimo Rainbow technology for SpO₂, SpCO and SpMet measurements is available as an option.

The carboxyhaemoglobin (SpCO) enables the non-invasive and immediate detection of heightened levels of carbon monoxide in the blood. This ensures rapid diagnosis and treatment of patients with carbon monoxide poisoning.

The SpMet measurement enables the non-invasive and immediate detection of high levels of methemoglobin in the blood.

NON-INVASIVE BLOOD PRESSURE MEASUREMENT

Arterial blood pressure measurement in children and adults, as well as a venous block function

INVASIVE BLOOD PRESSURE MEASUREMENT

Arterial or venous invasive blood pressure measurement with various pressure sensors

TEMPERATURE

Continuous measurement of the temperature through skin sensors or endo-cavity probes.

DIRECT DATA TRANSMISSION



INTERVENTION SITE

Recording of a 12-channel ECG on site and immediate transmission to the referring physician.




MEDICAL REGULATION

Reading and analysis of the data transmitted from the intervention site.



EMERGENCY TRANSPORT

The emergency doctor can implement the instructions immediately and thus help saving essential minutes.

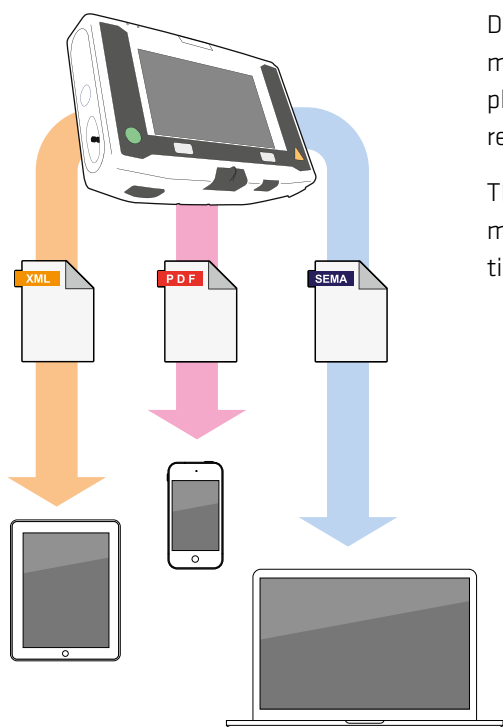

Direct data transmission ensures optimal care of the patient. The rescuer or physician can remain in contact with the regulating doctor.

The DGT7 / PGT7 enables data transmission during and after the intervention.

DURING THE INTERVENTION

- ❖ 12-lead ECG
- ❖ Trends
- ❖ Screenshots
- ❖ PDF reports
- ❖ Connection with ePCR

AFTER THE INTERVENTION

- ❖ Comprehensive intervention data including ECG and SpO₂ curves, events and trends.



CONNECTED AT ALL TIMES

- ❖ This device is compatible with the SCHILLER Update Server, which makes the updating through the network easier than ever before.
- ❖ Manage the entire configuration of the DGT7 / PGT7: the devices can be organised into groups, build configuration templates and be deployed easily to many devices at once.
- ❖ Device self-test results are automatically sent by e-mail, for a better pool management and immediate action in case of error.

- ❖ Thanks to its Bluetooth interface, the DGT7 / PGT7 offers the possibility to print Resting ECGs, live ECGs sequences, etc. wirelessly on a thermal printer.



DEFIGARD TOUCH 7



ARGUS PRO LIFEPOINT
The reusable CPR-Feedback-Sensor is easy to position.



ASSISTANCE IN DELIVERING CARDIAC MASSAGE

Thanks to the LifePoint CPR feedback sensor, cardiac massage monitoring is made easy. The following elements are displayed with various colors. Available information:

- Compression depth and rate
- Pressure relief at the end of each compression¹



MANUAL OR SEMI-AUTOMATIC DEFIBRILLATOR WITH UP TO 200 JOULES

The DEFIGARD Touch 7 can be operated in semi-automatic or manual mode (synchronously or asynchronously) with Multipulse Biowave (biphasic), which is achieved with exceptional efficacy and safety, as well as with extremely low energy input.

This makes this device ideal for all patient types. It identifies paediatric electrodes and adjusts the emitted energy automatically.



EXTERNAL PACEMAKER

- Transcutaneous stimulation via pads in Fix and Demand mode
- Wide range of settings

¹ Not valid for Great Britain, Germany and France

CONTAINS ALL FUNCTIONS



AED

MANUAL DEFIBRILLATOR

PACEMAKER

CPR FEEDBACK

COMPACT

LIGHTWEIGHT

TOUCH SCREEN

12 LEADS ECG

SPO₂, SPCO, SPMET

NIBP

ETCO₂ MAINSTREAM

ETCO₂ SIDESTREAM

IBP

TEMPERATURE

AIRCRAFT COMPATIBLE

AMBULANCE COMPATIBLE

WIFI

GSM

BLUETOOTH

USB

HIGH CONTRAST DISPLAY

NUMERIC VALUE VIEW



Americas
SCHILLER Americas Inc.
Doral, Florida 33172

North America:
Phone +1 786 845 06 20
Fax +1 786 845 06 02
sales@schilleramerica.com
www.schilleramerica.com

Latin America & Caribbean
Phone +1 305 591 11 21
Fax +1 786 845 06 02
sales@schillerla.ch
www.schillerla.ch



Asia
SCHILLER Asia-Pacific / Malaysia
52200 Kuala Lumpur, Malaysia
Phone +603 6272 3033
Fax +603 6272 2030
sales@schiller.com.my
www.schiller-asia.com



Austria
SCHILLER Handelsgesellschaft m.b.H.
A-4040 Linz
Phone +43 732 70 99
Fax +43 732 757 000
sales@schiller.at
www.schiller.at



China
Alfred Schiller (Beijing) Medical Equipment Co., Ltd.
100015 Beijing, China
Phone +86 010 52007020
Fax +86 010 52007020-8016
info@schillerchina.com
www.schillermedical.cn



Croatia
Schiller medicinski instrumenti d.o.o.
10000 Zagreb
Phone +385 1 309 66 59
Fax +385 1 309 66 60
info@schillerzg.hr
www.schiller.hr



France
SCHILLER Médical S.A.S.
F-67162 Wissembourg
Phone +33 3 88 63 36 00
Fax +33 3 88 63 36 49
info@schiller.fr
www.schiller-medical.com



France (distribution France)
SCHILLER France S.A.S.
F-77608 Bussy St Georges
Phone +33 1 64 66 50 00
Fax +33 1 64 66 50 10
contact@schillerfrance.fr
www.schiller-france.com



Germany
SCHILLER Medizintechnik GmbH
D-85622 Feldkirchen b. München
Phone +49 89 62 99 81 0
Fax +49 89 609 50 90
info@schillermed.de
www.schillermed.de



Hungary
SCHILLER Diamed Ltd.
H-1141 Budapest
Phone +36 1 383 4780 / +36 1 460 9491
Fax +36 1 383 4778
info@schillerhungary.hu
www.schillerhungary.hu



India
SCHILLER Healthcare India Pvt. Ltd.
Mumbai - 400 059, India
Phone +91 22 6152 3333 / 2920 9141
Fax +91 22 2920 9142
sales@schillerindia.com
www.schillerindia.com



Poland
SCHILLER Poland Sp. z o.o.
PL-02-729 Warszawa
Phone +48 22 843 20 89 / +48 22 647 35 90
Fax +48 22 843 20 89
schiller@schiller.pl
www.schiller.pl



Russia & C.I.S.
AO SCHILLER.RU
125124 Moscow, Russia
Phone +7 (495) 970 11 33
Fax +7 (495) 956 29 10
mail@schiller.ru
www.schiller.ru
www.schiller-cis.com



Serbia
SCHILLER d.o.o.
11010 Beograd
Phone +381 11 39 79 508
Fax +381 11 39 79 518
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
2310 Slovenska Bistrica
Phone +386 2 843 00 56
Fax +386 2 843 00 57
info@schiller.si
www.schiller.si



Spain
SCHILLER ESPAÑA, S.A.
E5-28232 Las Rozas/Madrid
Phone +34 91 713 01 76
Fax +34 91 355 79 33
schiller@schiller.es
www.schiller.es



Switzerland
SCHILLER-Reomed AG
CH-8953 Dietikon
Phone +41 44 744 30 00
Fax +41 44 740 37 10
sales@schiller-reomed.ch
www.schiller-reomed.ch



Turkey
SCHILLER TÜRKİYE
Okmeydanı-Sisli - İstanbul
Phone +90 212 210 8681 (pbx)
Fax +90 212 210 8684
info@schiller.com.tr
www.schiller-turkiye.com



United Kingdom
SCHILLER UK
Bellshill, ML4 3PR
Phone +44 1698 744 505
Fax +44 1698 744 474
sales@schilleruk.com
www.schilleruk.com



SCHILLER
The Art of Diagnostics