

medilog®AR


Efficient and reliable state-of-the-art Holter



SCHILLER
The Art of Diagnostics

medilog®AR

- ∴ Robust, shock- and splash-proof as well as easy to clean
- ∴ Flexible dual-battery concept for more than 14 days recording duration
- ∴ Zero seconds atrial fibrillation/atrial flutter detection based on true P-wave analysis

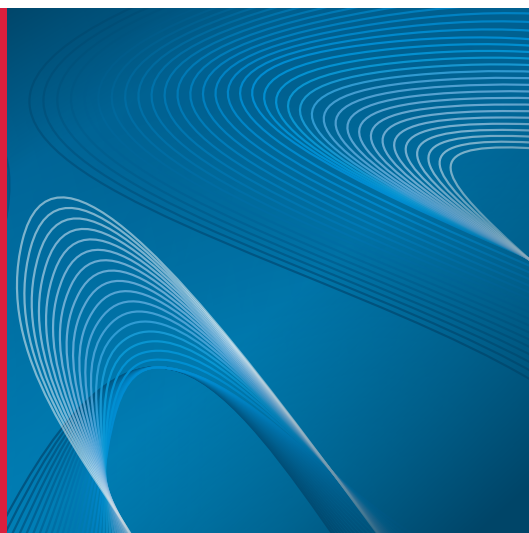

IDEALLY SUITED FOR CLINICAL ROUTINE

The shock- and splash-proof design makes the medilog®AR Holter highly suitable for the daily clinical routine. The Holter and its carrying case are robust and easy to clean.

TWO WEEKS OF POWER

Thanks to the dual-battery concept, the patients can be screened for more than 14 days without having to see the doctor to change the batteries.



Eco-friendly dual-battery concept (illustration)

ECO-FRIENDLY BATTERY:

By adjusting the signal resolution, the internal battery can run for 24/72 hours or more than 7 days. The high-quality internal battery is recharged via USB.

AAA SINGLE USE BATTERY:

For fastest turnaround times in high throughput clinics, use easily available standard AAA batteries.

MORE THAN 14 DAYS:

For maximum recording duration, combine a fully charged internal battery with an AAA battery to get more than 14 days with automated AF detection.

VARIOUS CARRYING OPTIONS

medilog®AR is so flat that it can easily be carried in a trouser pocket. With the carrying case, it can be worn on a belt or with a neck belt directly on the body.

MADE TO SUPPRESS ARTIFACTS

The medilog®AR Holter recorder features a superior resting ECG grade 32,000Hz sampling rate on 3 channels. This is used for true P-wave analysis, sophisticated artifact suppression and motion detection, enabling a speedy evaluation.

ZERO SECONDS ATRIAL FIBRILLATION/AFL DETECTION:

The medilog®AR scans for atrial fibrillation and atrial flutter by directly looking for P-waves, just as a cardiologist would.



RESPIRATION EPISODES:

Thanks to ECG-derived respiration recording, the medilog®AR is able to screen for potential respiratory episodes during sleep. The optional SpO₂ sensor, connected via Bluetooth, allows for additional respiratory information.



SpO₂ sensor, attached for the night

MUCH MORE THAN ARRHYTHMIA SCREENING

Locating Atrial Fibrillation in zero seconds, analysis of respiratory episodes, and high-resolution HRV are just some of the exceptional features of the medilog®AR. The analysis is carried out using the renowned medilog®DARWIN2 software.



Americas
SCHILLER Americas Inc.
 Doral, Florida 33172

North America:
 Phone +1 786 845 06 20
 Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com

Latin America & Caribbean
 Phone +1 305 591 11 21
 Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com



Asia
SCHILLER Asia-Pacific / Malaysia
 52200 Kuala Lumpur, Malaysia
 Phone +603 6272 3033
 Fax +603 6272 2030
sales@schiller.com.my
www.schiller-asia.com



Austria
SCHILLER Handelsgesellschaft m.b.H.
 A-4040 Linz
 Phone +43 732 70 90
 Fax +43 732 757 000
sales@schiller.at
www.schiller.at



China
Alfred Schiller (Beijing) Medical Equipment Co., Ltd.
 100015 Beijing, China
 Phone +86 010 52007020
 Fax +86 010 52007020-8016
info@schillerchina.com
www.schillermedical.cn



Croatia
SCHILLER d.o.o.
 10000 Zagreb
 Phone +385 1 309 66 59
 Fax +385 1 309 66 60
info@schillerzg.hr
www.schiller.hr



France
SCHILLER Médical
 F-67160 Wissembourg
 Phone +33 3 88 63 36 00
 Fax +33 3 88 63 36 49
info@schiller.fr
www.schiller-medical.com



France (distribution France)
SCHILLER France S.A.S.
 F-77608 Bussy St Georges
 Phone +33 1 64 66 50 00
 Fax +33 1 64 66 50 10
contact@schillerfrance.fr
www.schiller-france.com



Germany
SCHILLER Medizintechnik GmbH
 D-85622 Feldkirchen b. München
 Phone +49 89 62 99 81 0
 Fax +49 89 609 50 90
info@schillersed.de
www.schillersed.de



Hungary
SCHILLER Diamed Ltd.
 H-1141 Budapest
 Phone +36 1 383 4780 / +36 1 460 9491
 Fax +36 1 383 4778
info@schillerhungary.hu
www.schillerhungary.hu



India
SCHILLER Healthcare India Pvt. Ltd.
 Mumbai - 400 059, India
 Phone +91 22 6152 3333 / 2920 9141
 Fax +91 22 2920 9142
sales@schillerindia.com
www.schillerindia.com



Poland
SCHILLER Poland Sp. z o.o.
 PL-02-729 Warszawa
 Phone +48 22 843 20 89 / +48 22 647 35 90
 Fax +48 22 843 20 89
schiller@schiller.pl
www.schiller.pl



Russia & C.I.S.
AO SCHILLER.RU
 125040 Moscow, Russia
 Phone +7 (495) 970 11 33
 Fax +7 (495) 956 29 10
mail@schiller.ru
www.schiller.ru



Serbia
SCHILLER d.o.o.
 11010 Beograd
 Phone +381 11 39 79 508
 Fax +381 11 39 79 518
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
 2310 Slovenska Bistrica
 Phone +386 2 843 00 56
 Fax +386 2 843 00 57
info@schiller.si
www.schiller.si



Spain
SCHILLER ESPAÑA, S.A.
 ES-28232 Las Rozas/Madrid
 Phone +34 91 713 01 76
 Fax +34 91 355 79 33
schiller@schiller.es
www.schiller.es



Switzerland
SCHILLER-Reomed AG
 CH-8953 Dietikon
 Phone +41 44 744 30 00
 Fax +41 44 740 37 10
sales@schiller-reomed.ch
www.schiller-reomed.ch



Turkey
SCHILLER TÜRKİYE
 Okmeydani-Sisli - Istanbul
 Phone +90 212 210 8681 (pbx)
 Fax +90 212 210 8684
info@schiller.com.tr
www.schiller-turkiye.com



United Kingdom
SCHILLER UK
 Bellshill, ML4 3PR
 Phone +44 1698 744 505
 Fax +44 1698 744 474
sales@schilleruk.com
www.schilleruk.com

Device availability in your market is subject to regulatory approval.

