

Experience the freedom  
One ventilator, one circuit, multiple sites of care

LTV® 1200 ventilator (patients 5 kg to adult)

## The ventilation tools you need, throughout the continuum of care, for patients 5 kg and above.

At a fraction of the size of comparably equipped systems, the LTV® 1200 ventilator provides the ventilation tools you need, where you need them—critical care, emergency departments and emergency preparedness, patient transport, and long term care.

### A wide range of ventilation therapies for patients

The LTV® ventilator offers both invasive and noninvasive modes of ventilation in pressure control, Pressure Support, Volume Control and Spontaneous breath types. Studies show the use of noninvasive ventilation, when clinically appropriate, reduces incidents of VAP, improves weaning success and reduces total days of ventilation.

### Patient presets for quick initiation of ventilation


Choose infant, pediatric or adult to quickly configure initial ventilation, before optimizing settings for each patient's needs. This feature is proven as a simple and effective feature for emergency set-up.

### Unique Spontaneous Breathing Trial (SBT)

The SBT feature enhances the clinician's effectiveness during weaning with the potential to save time and money in the process. The SBT works using Rapid Shallow Breathing Index (RSBI) criteria to assess a patient's ability to be weaned. In addition, clinicians can optimize trial settings to each individual patient, ensuring the most favorable levels of support throughout the weaning process. The SBT feature allows the patient to synchronize their breathing pattern with pressure support and/or CPAP with variable rise time, flow termination and time termination.

### Oxygen resource management tool

The O<sub>2</sub> Conserve feature reduces oxygen consumption when the time and place of care require it. This feature extends the sometimes scarce resource of oxygen during transport or emergencies, and may be turned off when O<sub>2</sub> resources are again in full supply.

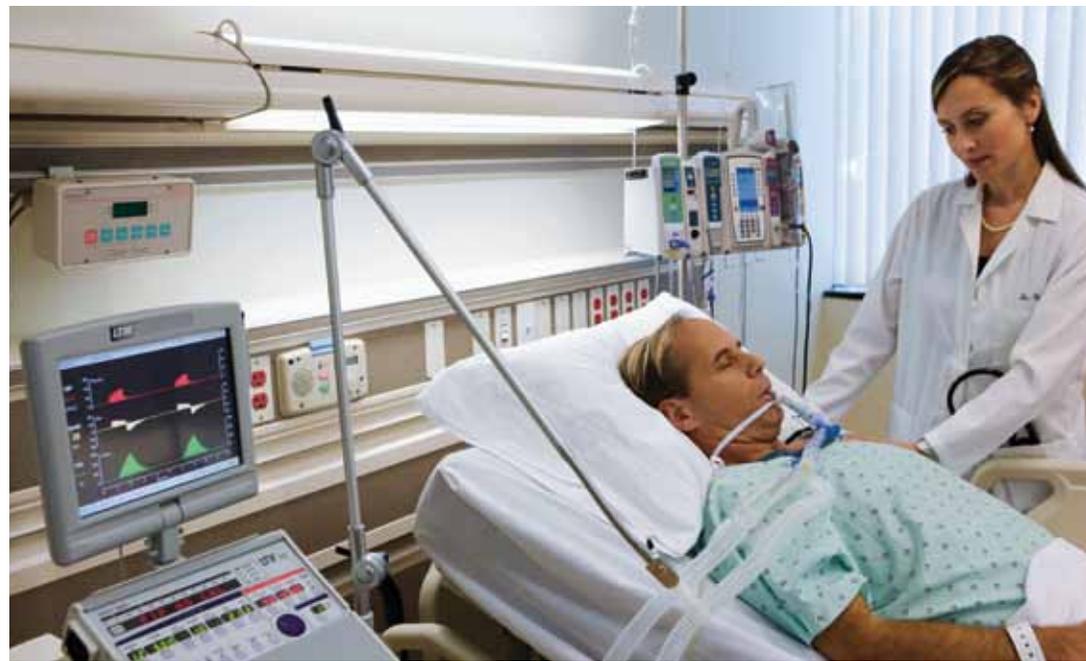

## Enhanced monitoring capabilities

Upgrade the system with the LTM™ Graphics Monitor, and the LTV® 1200 ventilator provides clinicians with even more comprehensive monitoring of critical patient-ventilator interaction. The graphics monitoring package displays real-time Pressure, Flow and Volume curves with adjustable cursors for accurate breath data measurement. The LTM™ also displays Flow/Volume loops and Volume/Pressure loops to enhance patient assessment. Patient data can be trended for up to 24 hours for clinical evaluation.

## Customizable system configurations

The LTV® 1200 ventilator system can be custom configured with a wide range of accessories to meet the needs of the most demanding care environments.

Accessories	Part Number
LTM™ Graphics Monitor package	17650-001
Floor stand	10611
Short crossbar	11451
Long crossbar	11452
Circuit support arm	11501
Humidifier pole	11453
Utility basket	11450
SprintPack lithium ion power system	19222-001
SprintPack bracket	19096-001
Transport pack for SprintPack	19102-001

## Specifications

Physical specifications with boots (approx.)	
Weight	14.5 lbs (6.5 kg)
Height	3.25 inches (8.4 cm)
Width	10.5 inches (27 cm)
Depth	13.5 inches (38 cm)
Variable alarms	
Apnea interval	10 to 60 seconds
High pressure limit	5 to 100 cmH <sub>2</sub> O
Low peak pressure	Off, 1 to 60 cmH <sub>2</sub> O
Low minute volume	Off, 0.1 to 99 liters
High PEEP	Off, 3 to 20 cmH <sub>2</sub> O
Low PEEP	Off, -3 to -20 cmH <sub>2</sub> O
High rate	Off, 5 to 80 bpm
Alarm volume	5 to 80 dBA at 1 meter

Controls	
Power	On/standby
Ventilation modes	Control, Assist/Control, SIMV/CPAP and NPPV
Breath types	VC, PC, PS and Spontaneous
Breath rate	0 to 80 bpm
Tidal volume	50 to 2000 ml
Inspiratory time	0.3 to 9.9 seconds (100 lpm)
PC/PS/Spont. flow	160 lpm
Pressure control	1 to 99 cmH <sub>2</sub> O
Pressure support	1 to 60 cmH <sub>2</sub> O
Sensitivity	Off, 1 to 9 lpm
% O <sub>2</sub>	21% to 100%
O <sub>2</sub> Flush	1 to 3 minutes
Low pressure O <sub>2</sub>	On/Off
PEEP/CPAP	0 to 20 cmH <sub>2</sub> O
Insp/Exp hold	6 second maximum
Manual breath	1 x current settings
Control lock	Easy or hard unlock options
Monitors and indicators	
Peak inspiratory pressure	0 to 120 cmH <sub>2</sub> O
Mean airway pressure	0 to 99 cmH <sub>2</sub> O
PEEP	0 to 99 cmH <sub>2</sub> O
Breath rate	0 to 250 bpm
Airway pressure display	-10 to 108 cmH <sub>2</sub> O
Exhaled tidal volume	0 to 4000 ml
Exhaled minute volume	0 to 99.9 liters
I:E Ratio	99:1 to 1:99
Calculated peak flow	10 to 100 lpm
AutoPEEP	0 to 100 cmH <sub>2</sub> O
Static compliance	1 to 999 ml/cmH <sub>2</sub> O
Patient effort	Green LED
Fixed alarms	
Disconnect/sense line	
External power low and lost	
High and low oxygen inlet pressure	
Internal battery low & empty	
Ventilator inoperative	

Extended features	
Spontaneous breathing trial (SBT)	
Ventilator presets (Infant, Pediatric and Adult)	
Variable rise time	
Variable flow termination	
Variable time termination	
Pressure control flow termination	
Leak compensation	
O <sub>2</sub> Conserve	
O <sub>2</sub> Cylinder duration calculation	
O <sub>2</sub> Flush	
Apnea backup ventilation	
Pneumatic specifications	
Pressure oxygen	40 to 80 PSIG (2.8 to 5.5 bar)
Low pressure source	< 80 lpm, < 10 PSIG
Power indicators	
External power	Green and amber LEDs
Battery charge status	Green, amber and red LEDs
Battery level	Green, amber and red LEDs
Power specifications	
Ventilator input 11 to 15 VDC	
AC Adaptor input 90 to 250 VAC	
50 to 60 Hz	
Environmental specifications	
Operating temperature	5° to 40° C (40 to 104 F)
Storage temperature	-20° to 60° C (-4 to 140 F)
Operating humidity	15% to 95% relative
Storage humidity	15% to 95% relative
Shock compliance	MIL-STD-810E
Vibration compliance	MIL-STD-810E
Standards and regulatory compliance	
cETLus Listed	
ISO 13485 registered (FM 75092)	
IEC 601-1-2	

For more information, please contact Customer Service within the USA at: 800.754.1914, and International: +1 763.398.8300.

LTV® 1200 product specifications are subject to change. Please consult your Sales Representative or CareFusion if you have any questions regarding updates to the product specifications.

CE 0086



Manufactured for:  
CareFusion  
22745 Savi Ranch Parkway  
Yorba Linda, CA 92887-4645 USA

EC REP

Authorized Representative:  
CareFusion Germany 234 GmbH  
Leibnizstrasse 7  
D-97204 Hoechberg, Germany

© 2010 CareFusion Corporation or one of its subsidiaries. All rights reserved. LTM and LTV are trademarks or registered trademarks of CareFusion Corporation or one of its subsidiaries. RC1328 (0610/4000)

18256-001 Rev. B

CareFusion  
Yorba Linda, CA

[carefusion.com](http://carefusion.com)

